


XVIVO Lung Cannula Set™

For connection of isolated donor lungs to an extracorporeal perfusion system

XVIVO Lung Cannula Set™ is a sterile, single-use set to be used to connect isolated donor lungs to an extracorporeal perfusion system for ex-vivo assessment. The cannulas create a tight connection, to allow an interface and provide circulatory support to the lung.

Features

- Built in pressure lines for reliable readings to avoid excessive hydrostatic pressure
- The left atrium cannula splints the left atrium open to create a reliable and consistent outflow drainage
- Easy to connect – with luer lock connections
- Colour coded to indicate PA (yellow) and LA (green) usage
- Designed to fit with XVIVO Organ Chamber™
- Material selected for optimal stability yet soft enough for running sutures
- Made of PVC and medical grade TPE. Not made with DEHP and latex.

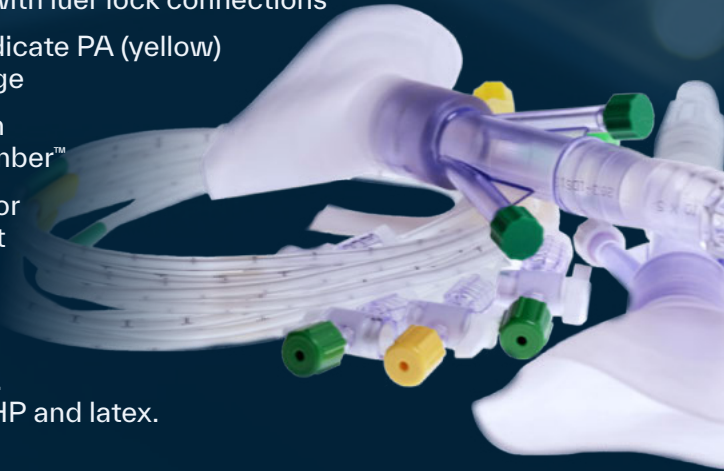

For US: The XVIVO Lung Cannula Set™ is only for use with XPS™ and is included in the XPS Disposable Lung Kit.

Order information:

Reference no:

19022

Product:

1 x XVIVO Lung Cannula Set™

XVIVO Perfusion AB, Mässans gata 10, SE-412 51 Göteborg, Sweden

Tel +46 31 788 21 50 | xvivogroup.com

XVIVO Lung Cannula Set™ is a trademark of XVIVO Perfusion AB.

© 2021 XVIVO Perfusion AB. All rights reserved.

Issued: 2021-10 Doc ID: MAR-16360-v.1.0

XVIVO